

INFLUENCE OF PHYSICAL AND CHEMICAL PARAMETERS ON EOSIN SPECTRAL DETERMINATIONS

Liljana Kola

University of Tirana, Faculty of Natural Sciences,
Center of Applied Nuclear Physics, P.O. Box 85, Tirana, Albania
kolaliljana@yahoo.it

The fluorescence ability of eosin enables its use as an artificial tracer in the water system studies. The problem deals with the application of eosin as a tracer to determine water movements within the karstic system and underground waters. The fluorescence intensity of fluorescent dyes in water samples depends on their physical and chemical properties, such as pH, temperature, and presence of oxidants, etc. In addition, the UV radiation may induce photochemical decomposition of eosin which can cause large measurement errors. This paper presents the results of the influence of these factors on eosin fluorescence intensity using the concentration and synchron scan methods. The method we have elaborated for this purpose made it possible to optimize procedures we use to analyze water samples for the presence of eosin and measure its content, even in trace levels by the means of the luminescent spectrometry

Key words: spectral determination; eosin; fluorescence intensity (I_F); synchron scan; artificial tracer

ВЛИЈАНИЕ НА ФИЗИЧКИТЕ И ХЕМИСКИТЕ ПАРАМЕТРИ ВРЗ СПЕКТРАЛНОТО ОПРЕДЕЛУВАЊЕ НА ЕОЗИН

Флуоресцентната способност на еозинот овозможува негова употреба како вештачки обележувач во студиите за водните системи. Ова својство се користи за следење и одредување на движењето на водата низ карсните предели и низ подземните текови. Интензитетот на флуоресценција на флуоресцентните бои во примероци на вода зависи од нивни физички и хемиски својства: рН, температура и присуство на оксиданси.. Освен тоа, УВ зрачењето може да индуцира фотохемиско распаѓање на еозинот, што може да доведе до големи грешки при негово мерење. Во овој труд се презентирани резултатите добиени со проучување на влијанието на споменатите фактори врз флуоресцентниот интензитет на еозин со користење на концентрациски и синхотронски скенирачки методи. Постапката која е образложена овозможува да се оптимизираат користените постапки за анализа на примероци вода во присуство на еозин. Содржината на еозин може да се определи дури и при негово присуство во траги, како што овде е направено со користење на луминисцентна спектрометрија.

Клучни зборови: спектрално определување; еозин; флуоресцентен интензитет (I_F); синхотронски скен; вештачки обележувач

1. INTRODUCTION

The management of aquatic systems requires a profound study them in order to assess their water usability as well as to find optimum and sustainable ways of water management [1, 2].

Modern and powerful techniques of tracers in hydrology are elaborated for this purpose. The in-

roduction of fluorescent and radioactive tracers has led to enormous methodological and instrumental developments [1, 2]. The problems they deal with when applied in studying various water systems might be grouped into different groups, according to where the groundwater is flowing to, from where it comes, whether underground hydraulic connections exist between different points

of the system, how the groundwater is flowing in(to) and/or through the system under study. [3].

The application of a spectrofluorometric technique in our Institute refers mainly to the study of water systems by using fluorescent dyes as artificial tracers. System's water is labelled according to some present plans and objectives and then their presence is monitored from time to time in various parts of the system under investigation. This information can be obtained by measuring the content of fluorescent dyes in direct, water samples or in active coal samples. In this paper we will present only the eosin results obtained from the experiments in water samples.

Eosin – $C_{20}H_6Br_4Na_2O_5$ is produced from the reaction of the fluoresceine and bromine [3, 4].

It is an anionic tracer not hazardous for humans, so it can be used in the study of water systems.

When measured in water samples eosin fluorescence intensity depends on different physical and chemical properties of the sample matrix – pH, temperature, oxidants' presence, etc. [3, 4, 5]. Their respective influence can be studied by luminescence spectrometry in order to determine optimal conditions for detecting the eosin maxima in trace levels.

2. EXPERIMENTAL

The fluorescence of water containing eosin was measured by means of Perkin-Elmer Spectral Fluorometer LS 55. We have previously analyzed some "blanks" sampled from different water systems, in order to assess the natural presence of the "eosin" fluorescence, the so called "background".

A special software package (FL WinLab) manages different application programs the instrument LS 55 offers. We have detected and measured tracer contents in water samples by using the Synchronous Scan and Concentration Applications. First the instrument validation should be realized through the Raman spectra (Raman Peak Wavelength, Raman Peak Intensity and Raman S/N ratio) from a sealed water cell [6] and a reference material (anthracene) [7]. Then the instrument stability is checked with regard to the fluorescence intensity by means of a reference material (anthracene) [5, 7]. Parameters were set in order to investigate the eosin content in water samples through synchronous and concentration applications. These

parameters were as follows: Excitation wavelength – 516 nm; Emission wavelength – 538 nm; Excitation slit – 10.0 nm; Emission slit – 10.0 nm; Delta lambda – 21 nm; Auto background subtract.

The characteristic peak wavelengths were verified beforehand by preliminary synchronous scans.

3. RESULTS AND DISCUSSION

Initially the instrument was calibrated with six eosin standard solutions by means of the calibration application for eosin determinations. The equation of the eosin calibration curve was:

$$I_F = 43.915 \cdot C + 0.019$$

$$R^2 = 0.9999$$

where:

- I – the fluorescence intensity
- C – eosin concentration (ppb)
- R^2 – the correlation coefficient

1. The influence of pH values on the eosin fluorescence

First, a measurement of the pH of the eosin standard solution and the respective fluorescence intensity in normal conditions was undertaken. Then the pH value was changed by adding HCl 0.05 mol/dm^{-3} through a micropipette, in order to create an acid medium, as well as EDTA-Na 0.05 mol/dm^{-3} , in order to create a basic medium. After each additional portion of the acid or basic solutions, the respective fluorescence of the sample was measured. The obtained results are shown in the following diagram (Fig. 1).

Figure 1 shows clearly the highest and the most stable values of the eosin fluorescence intensity are at pH values exceeding 5.58. By lowering pH below 5.50 the intensity declines rapidly. The eosin fluorescence intensity is considerably reduced when pH values of the sample are lowered below 5.58, reaching 5.84 % of its initial value for pH = 0.4. All this means in practice each time one is going to measure the water sample for the eosin presence can do it directly but in order to have correct measurements of the eosin content the sample matrix must be controlled before that pH > 5.6. The influence of pH on the eosin fluorescence intensity is a reversible process.

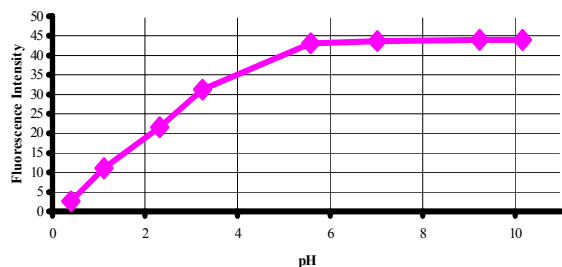


Fig. 1. Influence of pH on eosin fluorescence intensity

2. The influence of temperature on the eosin fluorescence intensity

When sample temperature is increased the excited molecules spend more energy in rotational and vibrational movements, which leads to heat production. It is expected therefore to measure lower fluorescence intensities [4, 5]. The temperature influence on fluorescence intensity of different dye tracers might be manifested differently. This phenomenon for the case of eosin has been studied.

First, the temperature and fluorescence intensity of the standard solution of 1 ppb eosin for normal conditions have been measured (the bold values in Table 1) and then we have changed temperature and measured respective fluorescence intensities. Table 1 shows the measurement results.

Table 1

Influence of temperature on intensity of eosin

No	Temperature °C	Intensity	Conc. ppb
1.	10.1	43.951	1.001
2.	12.9	43.461	0.989
3.	16.4	43.272	1.002
4.	20.5	43.968	1.001
5.	26.7	43.441	0.989
6.	29.8	43.756	0.996
7.	32.0	43.555	0.991
8.	37.9	43.671	0.994
9.	41.0	43.784	0.996

The measurement results demonstrate clearly that the eosin fluorescence is stable if temperature values of the water sample change from 10 to 41 °C.

3. The Influence of oxidants on eosin fluorescence intensity.

The fluorescence ability could be irreversibly destroyed following the action of oxidants, such as ozone, chloride, etc, on fluorescent substances.

Having this in view, we have focused our study on investigating the influence of the sodium hypochlorite on the eosin fluorescence intensity because it is usually used as a disinfectant in drinking water supply in our country. We have first measured the fluorescence intensity in the standard solution of 0.05 ppb eosin having no oxidizing agents. We have had equal amount of solution of 0.05 ppb eosin in eight flasks. Then we have added different portions of NaOCl in every flask: 100, 200, 300, 400, 500, 700, 900 μ l and measured respective values of the eosin fluorescence intensity. The obtained results demonstrate clearly the destructive effect of oxidizing agents on the eosin fluorescence in the water samples (Fig. 2).

Figure 2 shows plainly this effect. The first portion of the oxidizing substance added to an eosin sample decreased its fluorescence by 14.38 % of its initial intensity; 500 μ l of oxidizing agent in the same sample reduces fluorescence to 50 % of its initial value whereas 700 μ l of NaOCl reduces fluorescence to 76.45 % of its initial value. That means the oxidants destroy vigorously and irreversibly the fluorescent capacities of eosin.

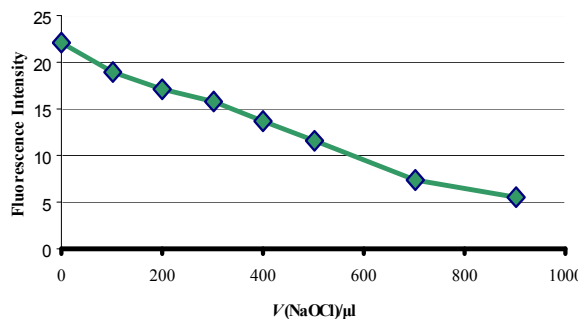


Fig. 2. Decay of eosin fluorescence in water solutions containing sodium hypochlorite

Further measurements carried out on this effect have shown that the fluorescence intensity is reduced in larger proportions when longer time periods of the oxidants presence in the samples are applied. The data in Table 2 are obtained by measuring the fluorescence intensity in different mo-

ments after having added sodium hypochlorite to the sample.

Table 2

The value of the I_F in different moments after having added NaOCl to the sample

Eosin st. sol. + NaOCl (ox.)	I_F $t = 0$ min	I_F $t = 30$ min	I_F $t = 60$ min	I_F $t = 90$ min
E + 0 μ l Ox	22.021	22.021	22.021	22.021
E + 100 μ l Ox.	18.862	12.832	8.576	4.906
E + 200 μ l Ox.	17.146	8.971	4.687	2.388
E + 300 μ l Ox.	15.755	4.922	1.795	0.676
E + 400 μ l Ox.	13.676	2.584	0.720	0.387
E + 500 μ l Ox.	11.547	1.360	0.552	0.264
E + 700 μ l Ox.	7.442	0.680	0.374	0.231
E + 900 μ l Ox.	5.612	0.530	0.252	0.135
E + 1000 μ l Ox.	5.118	0.473	0.244	0.135

In the flask containing eosin standard solution and 100 μ l NaOCl, only 68 % of the initial fluorescence has remained in it after 30 minutes while this number is lowered to 45.5 % after 60 minutes and to 26 % after 90 minutes. In the last flask containing eosin standard solution and 1000 μ l NaOCl, only 9.24 % of the initial fluorescence has remained in it after 30 minutes, while this value is lowered to 4.8 % after 60 minutes. After 90 minutes 97.4 % of the eosin initial fluorescence was destroyed.

That means the time factor should be kept always in view when working with water samples containing oxidants in studies using eosin as an artificial tracer.

4. The influence of the UV radiation on the eosin fluorescence intensity

The UV radiation may induce photochemical decomposition of eosin which can cause large measurement errors. To investigate this phenomenon, an eosin standard solution was prepared (conc. = 1 ppb, solvent: distilled water) and added in 2 colorless glass bottles and in 2 brown colored bottles. Then, they were put in two different situa-

tions so-called: 'in the light' and 'in the shadow'. One colorless glass bottle and one brown colored glass bottle were first put in the "in the light" and the other uncolored bottle and afterwards one brown colored one was put 'in the shadow'. The eosin fluorescence intensities of these four bottle were measured at the same time according to one schedule. At the beginning the measurements were frequent and after some times were less than before.

This process was followed for four months. The obtained results of the eosin photodecomposition in different conditions are shown in the following Tables and diagrams. Table 3 shows the decay of the fluorescence of eosin standard solution of 1 ppb in the colorless glass bottles exposed to daylight in the condition of so called 'in light'.

Table 3

Decay of the fluorescence of eosin aqueous solutions in colorless glass bottles, 'in light'

No	Exposition to daylight (day)	I_F of eosin	Eosin conc. (ppb)	Day of the measurement	Time of the measurement
1.	0	43.724	0.995	12. 03. 07	12: 00
2.	0.96	0.791	0.017	13. 03. 07	11: 30
3.	2.96	0.019	0.00	15. 03. 07	11:30
4	3.92	0.00	0.00	16. 03. 07	10:30

Table 3 shows how a 1 ppb eosin standard solution in colorless bottle decomposes when subjected to daylight. It is obvious from the results that the UV radiation decay immediately eosin. Light destroys 98.29 % of the total fluorescence intensity only after one day, whereas after two days exposure, the eosin fluorescence was destroyed 100 %.

Table 4 shows the results of the eosin fluorescence intensity (I_F) decay after exposure to daylight in the brown glass flask, 'in light'.

One can see clearly from the results of Table 4 that the fluorescence intensity of eosin after 3 days exposure to daylight in the brown glass flask, 'in light' is destroyed 91.53 % and completely after 42 days, so 39 days later than in the case of the colorless bottle in the same conditions. This proves that the brown flask delays the decomposition of eosin.

Table 4

Decay of the fluorescence intensity (I_F) of eosin in the brown glass bottle, 'in light'

No	Exposition to daylight (day)	I_F of eosin	Eosin conc. (ppb)
1.	0	43.724	0.995
2.	0.96	28.565	0.680
3.	2.96	3.704	0.084
4.	3.92	2.527	0.057
5.	6.92	0.245	0.005
6.	7.92	0.243	0.004
7.	8.92	0.195	0.004
8.	10.92	0.192	0.004
9.	13.92	0.187	0.004
10.	14.92	0.178	0.003
11.	15.92	0.157	0.003
12.	16.92	0.151	0.003
13.	17.92	0.086	0.002
14.	20.92	0.060	0.001
15.	37.92	0.030	0.001
16.	42.92	0.00	0.00
17.	50.92	0.00	0.00

The obtained results of the eosin fluorescence intensity (I_F) decay after exposure to daylight in colorless glass bottle, 'in shadow' are shown in Table 5.

Table 5

Decay of the eosin fluorescence intensity in the colorless glass bottle, 'in shadow'

No	I_F of eosin	Eosin conc. (ppb)	Exposition to daylight (day)
1.	43.724	0.995	0
2.	2.655	0.060	0.96
3.	0.234	0.005	2.96
4.	0.00	0.000	3.92

The results trace 5 show the decomposition of 1 ppb eosin standard solution in a colorless bottle when subjected to daylight. It is clear from the results that the UV radiation immediately destroys eosin in this case, too. Light destroys 93.93 % of the total fluorescence intensity only after one day.

The UV radiation destroyed 100 % of the eosin fluorescence after three days exposure, so one day later than in the case of the colorless bottle 'in light' (Tab. 3).

Table 6 shows the obtained results after exposure to daylight of 1 ppb eosin standard solution in a brown flask in condition 'in shadow'.

Table 6

Decay of the eosin fluorescence intensity in the brown glass bottle, 'in shadow'

No	I_F of eosin	Eosin conc (ppb)	Exposition to daylight (day)
1.	43.724	0.995	0
2.	41.013	0.933	0.96
3.	36.942	0.845	2.96
4.	35.456	0.807	3.92
5.	29.436	0.671	6.92
6.	29.072	0.662	7.92
7.	28.432	0.647	8.92
8.	26.056	0.593	10.92
9.	24.009	0.546	13.92
10.	23.346	0.531	14.92
11.	22.628	0.515	15.92
12.	22.051	0.502	16.92
13.	20.550	0.468	17.92
14.	17.385	0.395	20.92
15.	6.344	0.144	37.92
16.	6.242	0.142	42.92
17.	6.150	0.140	50.92
18.	6.066	0.138	59.92

The results of Table 6 show clearly that the decay of the eosin fluorescence intensity in the brown glass flask, 'in shadow' is more gradual than the decay in the other conditions. In this case, light destroys 86.3 % of the eosin total fluorescence intensity only after 60 days.

Figure 3 shows the dynamics of the eosin decomposition after exposure to daylight in four different conditions mentioned above.

It is obvious from the graphs that the brown glass and especially the brown flask delays the decomposition of eosin.

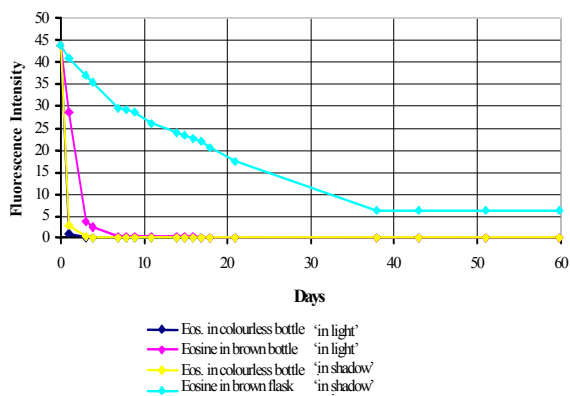


Fig. 3. Decay of the fluorescence of eosin in aqueous solutions exposed to different conditions of daylight

The results show clearly that eosin is extremely sensible to daylight. It can decay very quickly. That means one must take this influence into consideration when use eosin as an artificial tracer for water system studies. This study shows it should not be used in surface water tracing experiments but only in underground water studies.

4. CONCLUSIONS

- The maximal fluorescence of eosin in the aqueous solution is reached at the pH-value over 5.58. The fluorescence intensity of eosin decreases rapidly for $\text{pH} < 5.58$, but one can measure the eosin fluorescence intensity directly in water samples because their pH-values are usually over 6. However, the pH of water samples before eosin fluorescence measurements should be checked.

- One can observe that the eosin fluorescence intensity doesn't change for temperature range from 10 to 40 °C. That means special treatments should not be applied to the water samples and standards before eosin measurements.

- The oxidants destroy strongly and irreversibly the fluorescent capacities of eosin. This quenching effect on the eosin fluorescence is intensified at the concentration of oxidants water samples and exposure time of eosin to such agents.

For this reason, water must not be sampled from water supply network or it should be previously analyzed for probable presence of oxidants, and based on the results choose the sampling points. Otherwise the fluorescence measurements of such water samples could lead to erroneous results.

- The UV radiation induces a photochemical decomposition which leads to analytical errors and the destruction of eosin in the aqueous solution. That's why it is necessary to protect eosin standards and samples from the ambient light by storing them in brown glass and in the dark.

- Eosin is more affected by the ambient light than some other dyes and can decay very quickly. It should therefore not be used in surface water tracing experiments.

REFERENCES

- [1]. H. Behrens, *Tracer Methods in Isotope Hydrology*, IAEA, Vienna, 1983.
- [2]. H. Behrens, *Water Tracer Chemistry – Factor Determining Performance and Analytics of Tracers; 5th Int. Symp. on Underground Water Tracing*, Athens, 1986, pp. 121–133.
- [3]. W. Käss, *Tracing Technique in Geohydrology*, Berlin, 1994.
- [4]. R. Benischke, *Fluorescent Tracers in Hydrology*, Graz, 1991.
- [5]. L. Kola, *Optimization of Fluorometric Determinations in Water Samples and Active Coal Samples* (Albanian Printing), Ph.D. Thesis, Tirana University, 2008.
- [6]. Perkin Elmer, *Molecular Spectroscopy, LS 55 User's Guide*, Buckinghamshire, UK, 2000.
- [7]. Perkin Elmer, *Fluorescence Intensity Samples, C520-7440*, Booklet, 2000.